

ORION

ELECTRONIC PERSONAL SCALE

Model: ORION OS-06RD

USER'S GUIDE

ENGLISH

GB



Please read it carefully before use and save it

*Dear Customer,
Thank you for buying ORION products! Welcome to use!
Please read this User's Guide carefully before use and save it.
Best regards,
ORION ELECTRONICS Team
More details: <http://www.orion.ua>
Mail to us: support@orion.ua*

Care and Attention

- Place your scale on a hard flat surface will ensure the greatest accuracy and repeatability
- Always weigh yourself on the same scale placed on the same floor surface. Uneven floors may affect the readings.
- Always weigh yourself without clothing or footwear, and before meals. The first thing in the morning is the best time.
- Clean your scale with a damp cloth but do not allow water to get inside the scale. Do not use any chemical clearing agents.
- Do not allow your scale to become saturated with water as this can damage the electronics components.
- Do not store the scale in an upright position when not in use as this will drain the batteries.
- Remove the batteries if the scales are not going to be used for a long period of time.
- Treat your scale with care as it is a precision instrument, do not drop it or jump up or down on it.

Specifications

- Equipped with four high precision Strain Gauge sensors.
- Capacity 150kg/330lb/24st
- Division 100g/0.2lb/ 1/4st
- LCD Display
- Tap on to switch on.
- Automatic Zero resetting
- Automatic switch off
- Over load indicator
- Low battery indicator

Power Supply

- Lithium batteries are included (CR2032). Please remove the insulation sheet from the battery contacts before using. Replace the batteries when the low battery indicator shows.


Operating Instructions

- Step 1. Put the scale on a hard surface. (Avoiding carpets or any soft surface)
- Step 2. Switch on the scale with a slight tap on the glass platform near LCD display. The wait a few seconds the following will appear:



Note: The scale will not work properly if you step on the scale before it shows 0.0. If this happens please start from the beginning of step 2

- Step 3. Gently step on the scale. Put your feet evenly on the scale without moving and wait whilst the scale computes weight. Your correct weight will then be displayed in the following way: 80.3 Kg.

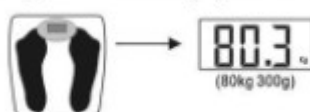

Automatic switch off

When you step off the scale, it will automatically switch off. (If no weighting is carried out while the zero is being displayed, the scale will automatically switch off after approximately 6 seconds).

Overload indicator

When the scale is overloaded (over 150kg/330lb/24st) the following will appear on the screen: ERR

Battery Change Indicator

When the batteries need replacing the following will appear on the screen: LO
Please open the back board of the scale, and then take the batteries carefully out of the battery compartment. Please push the metal sheet outwards by your finger nails or a sharp instrument.

The battery will pop up. Do not bend the contacts back while you do this, or they may break. Dispose of the old batteries in a safe place. Do not place on fire.

Subject to change without notice

Extra information on ORION company web-site www.orion.ua

We take your opinion and offers by E-mail support@orion.ua